QuantiChrom[™] α-L-Fucosidase Assay Kit (DFUC-100)

Quantitative Colorimetric Kinetic a-L-Fucosidase Activity Determination

DESCRIPTION

α-L-FUCOSIDASE (AFU) is an enzyme coded by the FUCA1 gene in humans and catalyzes the breakdown of L-Fucose. A genetic deficiency in this enzyme results in a neurovisceral storage disease, fucosidosis, which is characterized by the accumulation of fucose. Low serum activity of fucosidase has also been linked to ovarian caricinoma. Elevated fucosidase serum activity has been observed in patients with diabetes. hyperthyroidism, cirrhosis, and hepatitis. Increased activity has been associated with lung, breast, stomach, ovary, uterus, and liver carcinomas.

BioAssay Systems' non-radioactive, colorimetric AFU assay is based on the cleavage of 4-nitrophenol from the synthetic substrate. Nitrophenol becomes intensely colored after addition of the stop reagent. The increase in absorbance at 405 nm after addition of the stop reagent is directly proportional to the enzyme activity.

KEY FEATURES

High sensitivity and wide linear range. Linear detection range (10 µL sample): 1 to 100 U/L for a 20 minute reaction.

Homogeneous and simple procedure. Simple add-mix-read procedure allows reliable quantitation of fucosidase activity within 20 minutes.

Robust and amenable to HTS. All reagents are compatible with highthroughput liquid handling instruments. Can be readily automated to measure thousands of samples per day.

APPLICATIONS

α-L-Fucosidase activity determination in biological samples (e.g. plasma, serum, tissue, cell lysate, etc.)

KIT CONTENTS (100 TESTS IN 96-WELL PLATES)

Substrate Buffer: 10 mL Stop Reagent: 12 mL

1 mL 12.5 mM Nitrophenol Standard:

Storage conditions. The kit is shipped at room temperature. Store all components at -20°C upon receiving. Shelf life: 6 months after receipt.

Precautions: reagents are for research use only. Normal precautions for laboratory reagents should be exercised while using the reagents. Please refer to Material Safety Data Sheet for detailed information.

PROCEDURES

This assay is based on a kinetic reaction. To ensure identical incubation time, addition of Substrate and Stop Reagent to samples should be quick and mixing should be brief but thorough. Use of a multi-channel pipettor is recommended.

Sample Preparation: Serum and plasma can be assayed directly.

Tissue: Prior to dissection, rinse tissue in phosphate buffered saline (pH 7.4) to remove blood. Homogenize tissue (50 mg) in ~200 µL buffer containing 50 mM potassium phosphate (pH 7.5). Centrifuge at 10,000 × g for 15 min at 4°C. Remove supernatant for assay.

Cell Lysate: Collect cells by centrifugation at 2,000 × g for 5 min at 4°C. For adherent cells, do not harvest cells using proteolytic enzymes; rather use a rubber policeman. Homogenize or sonicate cells in an appropriate volume of cold buffer containing 50 mM potassium phosphate (pH 7.5). Centrifuge at 10,000 x g for 15 min at 4°C. Remove supernatant for assav.

All samples can be stored at -80 to -20°C for at least one month.

Reagent Preparation: Equilibrate Substrate Buffer to desired reaction temperature (e.g. 25°C or 37°C).

Standard Preparation: Mix 10 μL of 12.5 mM nitrophenol standard with 490 μL dH₂O to make 250 μM standard.

No	250 μM STD + dH2O	Vol (μL)	Nitrophenol (µM)
1	200 μL + 0 μL	200	250
2	120 μL + 80 μL	200	150
3	60 μL + 140 μL	200	75
4	0 μL + 200 μL	200	0

Reaction Preparation:

- 1. Transfer 100 μL of each standard (OD_{STD}) into wells of a clear flat bottom 96-well plate.
- 2. Transfer 20 μL of each sample into separate wells. Add 80 μL Substrate to each sample well. Tap plate briefly to mix.
- 3. Incubate at 25°C or desired temperature for 20 minutes. Add 100 µL of Stop Reagent to each well. Tap plate briefly to mix.

Note: If your sample is colored or opaque, then a sample blank (ODBLANK) will be needed. Add 20 μ L of sample to a well, and add 80 μ L of dH₂O. After incubation add 100 µL Stop Reagent.

CALCULATION

Subtract blank OD (water, #4) from the standard OD values and plot the ΔOD against standard concentrations. Determine the Slope and use the following equation to calculate α -Fucosidase activity:

AFU Activity =
$$\frac{\text{OD}_{\text{SAMPLE}} - \text{OD}_{\text{BLANK}}}{\text{Time • Slope}} \times \frac{\text{Reaction Vol } (\mu L)}{\text{Sample Vol } (\mu L)} \times n$$

$$= \frac{\text{OD}_{\text{SAMPLE}} - \text{OD}_{\text{BLANK}}}{\text{Slope}} \times \frac{1}{4} \times n \quad \text{(U/L)}$$


where OD_{SAMPLE} is the OD_{405nm} value for each sample and OD_{BLANK} is the OD_{405nm} value of the water (standard #4) or the sample blank if one was used. Slope is the slope of the linear regression fit of the standard points and Time is the reaction time (20 min). Reaction Vol and Sample Vol are 100 μ L and 20 μ L, respectively. *n* is the dilution factor.


Unit definition: 1 Unit (U) of AFU will catalyze the conversion of 1 µmole of 4-Nitrophenyl- α -L-fucopyranoside to 4-Nitrophenol and α -L-Fucose per min at 25°C and pH 5.3.

Note: If sample AFU activity exceeds 100 U/L, either use a shorter reaction time or dilute samples in water and repeat the assay. For samples with AFU activity < 5 U/L, the incubation time can be extended up to 40 minutes for greater sensitivity.

MATERIALS REQUIRED, BUT NOT PROVIDED

Pipetting devices and accessories (e.g. multi-channel pipettor), clear flatbottom 96-well plates (e.g. VWR cat# 82050-760), centrifuge tubes and plate reader.

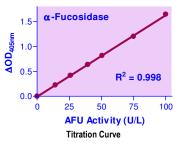

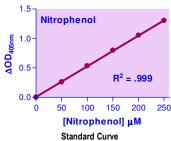

LITERATURE

- 1. Giardina M.G. et al. (1992) Serum alpha-L-fucosidase. A useful marker in the diagnosis of hepatocellular carcinoma. Cancer. 70:1044-1048
- 2. Fernandez-Rodriguez AD. Et al. (2000) Value of the serum alpha-Lfucosidase activity in the diagnosis of colorectal cancer. Oncology 59(4): 310-316.
- 3. Alhadeff JA et al. (1975) Human liver alpha-L-fucosidase. Purification, characterization, and immunochemical studies. Journal Biol. Chem. 250: 7106-7113.